

- L** LABORATORY
- P** PROCESS
- S** SOFTWARE
- A** AUTOMATION



**SCHMIDT
HAENSCH**
innovators by tradition since 1864


AQUISYS 2008

21 CFR Part 11 compliant Software




Aquisys 2008

Optional to our high precision measuring instruments, Schmidt + Haensch offers you the completely 21 CFR Part 11 compliant software **AQUISYS 2008**. This software enables you to control your polarimeters, refractometers and density meters with your PC.

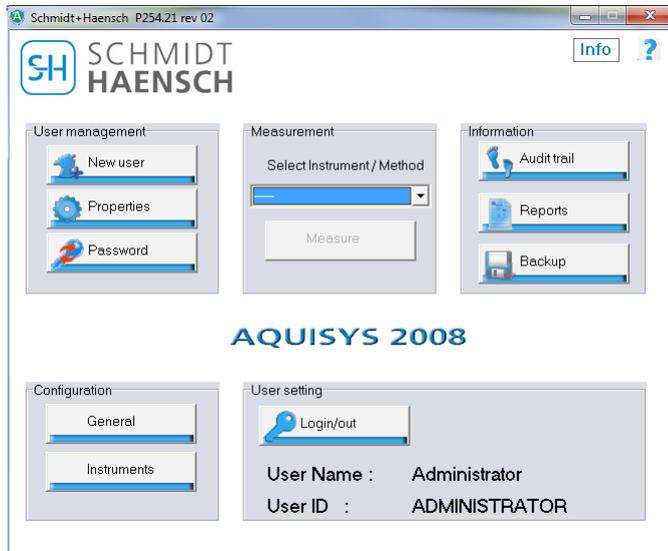


Aquisys 2008 Main Screen

- The integrated, though independent program "VALIDATE" checks the integrity of Aquisys before every run
- Remote operation of SCHMIDT+HAENSCH's laboratory instruments (Polartronic, Saccharomat and ATR) with your PC
- Audit trail function (history logging, one or all users divided by event type)
- Report generation
- Electronic signature which is unmistakably attributed to a person
- The system is able to generate accurate and complete paper prints of the electronic records
- "Manual signing" area at printouts
- Encrypted electronic data storage and archiving

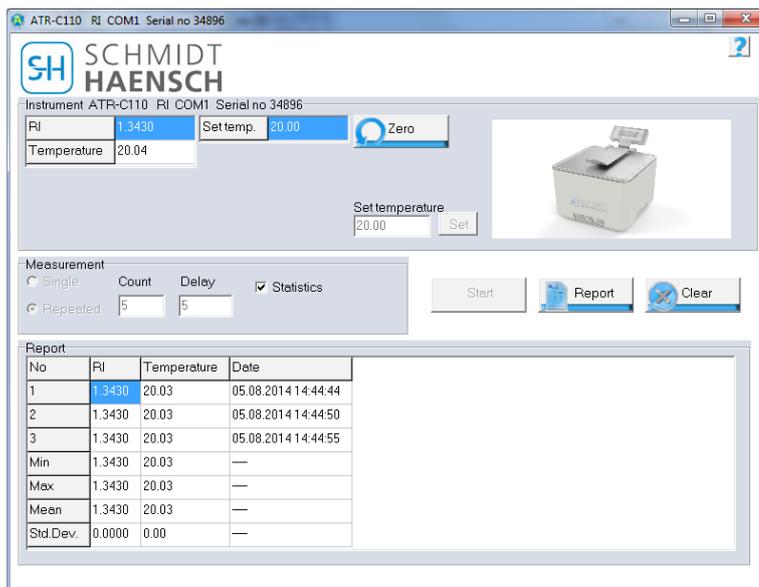

Aquisys 2008 Measuring Menu

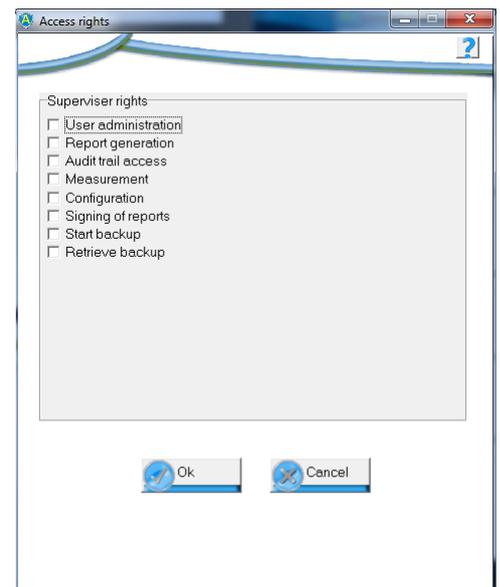

Aquisys 2008 Access Rights Menu

- Backup (of complete settings, audit trail and reports) function -> storing as .zip-files
- Reload backups
- User management with individually configurable access rights for every user
- Logging of electronic recordings and modifications including user information as well as date and time stamp
- Explicit identification of the connected measuring instruments
- Passwords with at least seven characters and one alphabetic character
- Encrypted remote mode
- Single and repeated measurements